



Solution Kit

Sinus Combination Kit

Introduction

- Composition of the essential tools for sinus elevation techniques
- An integrated sinus kit, considering user's convenience

Feature

- The sinus integration kit consists of tools which enable the sinus crestal approach technique & sinus lateral approach technique.




Direction for use


For Sinus Crestal Approach Kit




Initial Drill

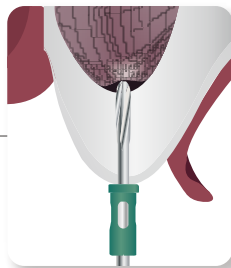

Sinus Reamer

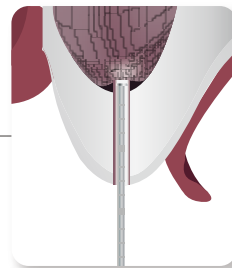



Depth Gauge

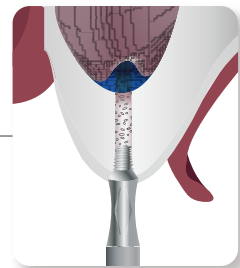


Aqua Lift Tip

For Sinus Lateral Approach Kit

Case1

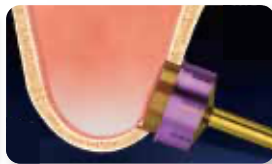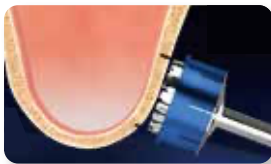

Case2

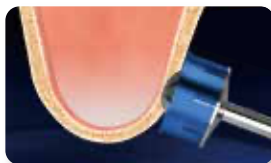

Sinus Combination Kit

Sinus Crestal Initial Drill



Ø1.8 Ø2.3

D	Ø1.8	Ø2.3
-	PD-18	PD-23

- Connecting with stopper, use for drilling until the cortical bone nearby sinus area.
 - Available to connect with stopper.
 - Same tool in the Sinus Crestal Approach Kit
- [Recommended R.P.M. : 800 ~ 1,000 R.P.M.]

Sinus Crestal Reamer



Ø2.8 Ø3.3 Ø3.8

D	Ø2.8	Ø3.3	Ø3.8
-	SRM-28	SRM-33	SRM-38

- Connecting with stopper, use for drilling where the initial drill is perforated until the cortical bone.
 - Even the reamer is reached to the sinus wall, the sinus wall isn't perforated.
 - Available to connect with stopper.
 - Same tool in the Sinus Combination Kit & Simple Sinus Crestal Reamer Kit
- Both kits have Ø 4.2 diameter of reamer (SRM-42) as well.
- [Recommended R.P.M. : 800 ~ 1,000 R.P.M.]

Sinus Crestal Diamond Bur



D	Ø2.8
-	DAR-28

- Even the bur is reached to the sinus wall, the sinus wall isn't perforated.
 - May be used for the final drill on the surgery on confirming sinus membrane drilling.
 - Available to connect with stopper.
 - Same tool in the Sinus Crestal Approach Kit & Simple Sinus Crestal Reamer Kit
- [Recommended R.P.M. : 1,000 ~ 1,200 R.P.M.]

Sinus Crestal Stopper



L	1.0mm	2.0mm	3.0mm	4.0mm	5.0mm	6.0mm	7.0mm	8.0mm
-	SCA-ST01	SCA-ST02	SCA-ST03	SCA-ST04	SCA-ST05	SCA-ST06	SCA-ST07	SCA-ST08

- Possible of safe drilling after connecting with drills or reamers.
- Same tool in the Sinus Crestal Approach Kit. But, in that kit, all ten kinds of stopper are composed.
- Same tool in the Simple Sinus Crestal Reamer Kit.
But, in that kit, five kinds of stopper (from 1.0mm to 5.0mm) are composed only.

Depth Gauge



-	DG-01
---	-------

- After drilling, used for the depth measurement or check whether the sinus wall is perforated.
- Available to connect with stopper.
- Same tool in the Sinus Crestal Approach Kit.

Aqua Tip



Ø2.8

Ø3.3

D	Ø2.8	Ø3.3
-	ALT-28	ALT-33

- Use for elevating sinus wall by water pressure
- After drilling by sinus crestal reamers, using water pressure elevate sinus wall safely.
- Aqua tip with Ø 2.8 Diameter is used when the sinus crestal reamers (Ø 2.8 ~ Ø 3.8) are perforated.
- Aqua tip with Ø 3.3 Diameter is used when the sinus crestal reamers (Ø 3.3 ~ Ø 4.2) are perforated.
- Same tool in the Sinus Crestal Approach Kit.

Silicone Tube



D	Ø3.0
-	AT-02

- Connecting with Aqua Tip, insert water into silicone tube to make a water pressure
- Same tool in the Sinus Crestal Approach Kit.

Sinus Lateral Core Drill



Ø7.0

D	Ø7.0
-	CD-70

- Connecting with stopper, use for perforating the lateral window
- Available to connect with stopper.
- Same tool in the Sinus Lateral Approach Kit & Simple Sinus Lateral Reamer Kit.
[Recommended R.P.M. : 1,000 ~ 1,200 R.P.M.]

Sinus Lateral Diamond Core Reamer



Ø7.0

D	Ø7.0
-	DCR-70

- Help to perforate safely after perforating the lateral window by core drill.
- Possible to use to check whether the lateral sinus wall is perforated.
- Available to connect with stopper.
- Same tool in the Sinus Lateral Approach Kit & Simple Sinus Lateral Reamer Kit
[Recommended R.P.M. : 1,200 ~ 1,500 R.P.M.]

Sinus Lateral Reamer



Ø8.0

Ø6.5/8.0

D	Ø8.0	Ø6.5/Ø8.0
-	SRM-80	SRM-80T

- Connecting with stopper, use for perforating the lateral window.
- The lateral reamer of Ø 6.5 / Ø 8.0 is used for flat perforation or window opening extension.
- Available to connect with stopper.
- Same tool in the Sinus Lateral Approach Kit.
But, in that Kits, Ø 6.5 diameter of reamer (SRM-65) is composed as well.
- Same tool in the Simple Sinus Lateral Reamer Kit
[Recommended R.P.M. : 1,000 ~ 1,200 R.P.M.]

Sinus Lateral Diamond Reamer



Ø6.5/8.0

D	Ø6.5 / Ø8.0
-	SDR-80T

- Help to perforate safely after perforating the lateral window by sinus lateral reamer.
- Possible to use to check whether the lateral sinus wall is perforated.
- Available to connect with stopper.
- Same tool in the Sinus Lateral Approach Kit & Simple Sinus Lateral Reamer Kit.
[Recommended R.P.M. : 1,200 ~ 1,500 R.P.M.]

Sinus Lateral Stopper



1.0mm



1.5mm



2.0mm



2.5mm



3.0mm

L	1.0mm	1.5mm	2.0mm	2.5mm	3.0mm
-	SCL-ST10	SCL-ST15	SCL-ST20	SCL-ST25	SCL-ST30

- Possible of safe perforation after connecting with drills or reamers
- Same tool in the Sinus Lateral Approach Kit & Simple Sinus Lateral Reamer Kit

Side Drill



Ø3.0

D	Ø3.0
-	SID-30

- Use for extension of sinus window, after perforation
- Same tool in the Sinus Lateral Approach Kit
[Recommended R.P.M. : 1,000 ~ 1,200 R.P.M.]

Sinus Curette



Nº	#1	#2	#3	#4	#5
-	7091	7092	7093	7094	7095

- After the sinus window is opened, use for elevating sinus wall manually.
- Same tool in the Sinus Lateral Approach Kit